# World Journal of *Clinical Cases*

World J Clin Cases 2020 October 6; 8(19): 4280-4687





Published by Baishideng Publishing Group Inc

W J C C World Journal of Clinical Cases

# Contents

# Semimonthly Volume 8 Number 19 October 6, 2020

# **OPINION REVIEW**

4280 Role of monoclonal antibody drugs in the treatment of COVID-19 Ucciferri C, Vecchiet J, Falasca K

# **MINIREVIEWS**

- 4286 Review of simulation model for education of point-of-care ultrasound using easy-to-make tools Shin KC, Ha YR, Lee SJ, Ahn JH
- 4303 Liver injury in COVID-19: A minireview Zhao JN. Fan Y. Wu SD

# **ORIGINAL ARTICLE**

# **Case Control Study**

4311 Transanal minimally invasive surgery vs endoscopic mucosal resection for rectal benign tumors and rectal carcinoids: A retrospective analysis

Shen JM, Zhao JY, Ye T, Gong LF, Wang HP, Chen WJ, Cai YK

4320 Impact of *mTOR* gene polymorphisms and gene-tea interaction on susceptibility to tuberculosis Wang M, Ma SJ, Wu XY, Zhang X, Abesig J, Xiao ZH, Huang X, Yan HP, Wang J, Chen MS, Tan HZ

# **Retrospective Cohort Study**

4331 Establishment and validation of a nomogram to predict the risk of ovarian metastasis in gastric cancer: Based on a large cohort

Li SQ, Zhang KC, Li JY, Liang WQ, Gao YH, Qiao Z, Xi HQ, Chen L

# **Retrospective Study**

4342 Predictive factors for early clinical response in community-onset Escherichia coli urinary tract infection and effects of initial antibiotic treatment on early clinical response

Kim YJ, Lee JM, Lee JH

- 4349 Managing acute appendicitis during the COVID-19 pandemic in Jiaxing, China Zhou Y, Cen LS
- 4360 Clinical application of combined detection of SARS-CoV-2-specific antibody and nucleic acid Meng QB, Peng JJ, Wei X, Yang JY, Li PC, Qu ZW, Xiong YF, Wu GJ, Hu ZM, Yu JC, Su W
- Prolonged prothrombin time at admission predicts poor clinical outcome in COVID-19 patients 4370 Wang L, He WB, Yu XM, Hu DL, Jiang H



World Journal of Clinica Contents					
Conter	Semimonthly Volume 8 Number 19 October 6, 2020				
4380	Percutaneous radiofrequency ablation is superior to hepatic resection in patients with small hepatocellular carcinoma				
	Zhang YH, Su B, Sun P, Li RM, Peng XC, Cai J				
4388	Clinical study on the surgical treatment of atypical Lisfranc joint complex injury				
	Li X, Jia LS, Li A, Xie X, Cui J, Li GL				
4400	Application of medial column classification in treatment of intra-articular calcaneal fractures				
	Zheng G, Xia F, Yang S, Cui J				
	Clinical Trials Study				
4410	Optimal hang time of enteral formula at standard room temperature and high temperature				
	Lakananurak N, Nalinthassanai N, Suansawang W, Panarat P				
	META-ANALYSIS				
4416	Meta-analysis reveals an association between acute pancreatitis and the risk of pancreatic cancer				
	Liu J, Wang Y, Yu Y				
	SCIENTOMETRICS				
4431	Global analysis of daily new COVID-19 cases reveals many static-phase countries including the United States potentially with unstoppable epidemic				
	Long C, Fu XM, Fu ZF				
	CASE REPORT				
4443	Left atrial appendage aneurysm: A case report				
	Belov DV, Moskalev VI, Garbuzenko DV, Arefyev NO				
4450	Twenty-year survival after iterative surgery for metastatic renal cell carcinoma: A case report and review of literature				
	De Raffele E, Mirarchi M, Casadei R, Ricci C, Brunocilla E, Minni F				
4466	Primary rhabdomyosarcoma: An extremely rare and aggressive variant of male breast cancer				
	Satală CB, Jung I, Bara TJ, Simu P, Simu I, Vlad M, Szodorai R, Gurzu S				
4475	Bladder stones in a closed diverticulum caused by Schistosoma mansoni: A case report				
	Alkhamees MA				
4481	Cutaneous ciliated cyst on the anterior neck in young women: A case report				
	Kim YH, Lee J				
4488	Extremely rare case of successful treatment of metastatic ovarian undifferentiated carcinoma with high dose combination cytotoxic chemotherapy: A case report				
	Kim HB Lee HI Hong R Park SG				

Kim HB, Lee HJ, Hong R, Park SG



	World Journal of Clinical Cases					
Conter	nts Semimonthly Volume 8 Number 19 October 6, 2020					
4494	Acute amnesia during pregnancy due to bilateral fornix infarction: A case report					
	Cho MJ, Shin DI, Han MK, Yum KS					
4499	Ascaris-mimicking common bile duct stone: A case report					
	Choi SY, Jo HE, Lee YN, Lee JE, Lee MH, Lim S, Yi BH					
4505	Eight-year follow-up of locally advanced lymphoepithelioma-like carcinoma at upper urinary tract: A case report					
	Yang CH, Weng WC, Lin YS, Huang LH, Lu CH, Hsu CY, Ou YC, Tung MC					
4512	Spontaneous resolution of idiopathic intestinal obstruction after pneumonia: A case report					
	Zhang BQ, Dai XY, Ye QY, Chang L, Wang ZW, Li XQ, Li YN					
4521	Successful pregnancy after protective hemodialysis for chronic kidney disease: A case report					
	Wang ML, He YD, Yang HX, Chen Q					
4527	Rapid remission of refractory synovitis, acne, pustulosis, hyperostosis, and osteitis syndrome in response to the Janus kinase inhibitor tofacitinib: A case report					
	Li B, Li GW, Xue L, Chen YY					
4535	Percutaneous fixation of neonatal humeral physeal fracture: A case report and review of the literature					
	Tan W, Wang FH, Yao JH, Wu WP, Li YB, Ji YL, Qian YP					
4544	Severe fundus lesions induced by ocular jellyfish stings: A case report					
	Zheng XY, Cheng DJ, Lian LH, Zhang RT, Yu XY					
4550	Application of ozonated water for treatment of gastro-thoracic fistula after comprehensive esophageal squamous cell carcinoma therapy: A case report					
	Wu DD, Hao KN, Chen XJ, Li XM, He XF					
4558	Germinomas of the basal ganglia and thalamus: Four case reports					
	Huang ZC, Dong Q, Song EP, Chen ZJ, Zhang JH, Hou B, Lu ZQ, Qin F					
4565	Gastrointestinal bleeding caused by jejunal angiosarcoma: A case report					
	Hui YY, Zhu LP, Yang B, Zhang ZY, Zhang YJ, Chen X, Wang BM					
4572	High expression of squamous cell carcinoma antigen in poorly differentiated adenocarcinoma of the stomach: A case report					
	Wang L, Huang L, Xi L, Zhang SC, Zhang JX					
4579	Therapy-related acute promyelocytic leukemia with FMS-like tyrosine kinase 3-internal tandem duplication mutation in solitary bone plasmacytoma: A case report					
	Hong LL, Sheng XF, Zhuang HF					
4588	Metastasis of esophageal squamous cell carcinoma to the thyroid gland with widespread nodal involvement: A case report					
	Zhang X, Gu X, Li JG, Hu XJ					

World Journal of Clinical Cases				
Conter	semimonthly Volume 8 Number 19 October 6, 2020			
4595	Severe hyperlipemia-induced pseudoerythrocytosis - Implication for misdiagnosis and blood transfusion: A case report and literature review			
	Zhao XC, Ju B, Wei N, Ding J, Meng FJ, Zhao HG			
4603	Novel brachytherapy drainage tube loaded with double 125I strands for hilar cholangiocarcinoma: A case report			
	Lei QY, Jiao DC, Han XW			
4609	Resorption of upwardly displaced lumbar disk herniation after nonsurgical treatment: A case report			
	Wang Y, Liao SC, Dai GG, Jiang L			
4615	Primary hepatic myelolipoma: A case report and review of the literature			
	Li KY, Wei AL, Li A			
4624	Endoscopic palliative resection of a giant 26-cm esophageal tumor: A case report			
	Li Y, Guo LJ, Ma YC, Ye LS, Hu B			
4633	Solitary hepatic lymphangioma mimicking liver malignancy: A case report and literature review			
	Long X, Zhang L, Cheng Q, Chen Q, Chen XP			
4644	Intraosseous venous malformation of the maxilla after enucleation of a hemophilic pseudotumor: A case report			
	Cai X, Yu JJ, Tian H, Shan ZF, Liu XY, Jia J			
4652	Intravesically instilled gemcitabine-induced lung injury in a patient with invasive urothelial carcinoma: A case report			
	Zhou XM, Wu C, Gu X			
4660	Bochdalek hernia masquerading as severe acute pancreatitis during the third trimester of pregnancy: A case report			
	Zou YZ, Yang JP, Zhou XJ, Li K, Li XM, Song CH			
4667	Localized primary gastric amyloidosis: Three case reports			
	Liu XM, Di LJ, Zhu JX, Wu XL, Li HP, Wu HC, Tuo BG			
4676	Displacement of peritoneal end of a shunt tube to pleural cavity: A case report			
	Liu J, Guo M			
4681	Parathyroid adenoma combined with a rib tumor as the primary disease: A case report <i>Han L, Zhu XF</i>			

# Contents

Semimonthly Volume 8 Number 19 October 6, 2020

# **ABOUT COVER**

Peer-reviewer of World Journal of Clinical Cases, Prof. Adrián Ángel Inchauspe, obtained his MD in 1986 from La Plata National University (Argentina), where he remained as Professor of Surgery. Study abroad, at the Aachen and Tubingen Universities in Germany in 1991, led to his certification in laparoscopic surgery, and at the Louis Pasteur University in Strasbourg France, led to his being awarded the Argentine National Invention Award in 1998 for his graduate work in tele-surgery. He currently serves as teacher in the Argentine Acupuncture Society, as Invited Foreigner Professor at the China National Academy of Sciences and Hainan Medical University, and as editorial member and reviewer for many internationally renowned journals. (L-Editor: Filipodia)

# AIMS AND SCOPE

The primary aim of World Journal of Clinical Cases (WJCC, World J Clin Cases) is to provide scholars and readers from various fields of clinical medicine with a platform to publish high-quality clinical research articles and communicate their research findings online.

WJCC mainly publishes articles reporting research results and findings obtained in the field of clinical medicine and covering a wide range of topics, including case control studies, retrospective cohort studies, retrospective studies, clinical trials studies, observational studies, prospective studies, randomized controlled trials, randomized clinical trials, systematic reviews, meta-analysis, and case reports.

# **INDEXING/ABSTRACTING**

The WJCC is now indexed in Science Citation Index Expanded (also known as SciSearch®), Journal Citation Reports/Science Edition, PubMed, and PubMed Central. The 2020 Edition of Journal Citation Reports® cites the 2019 impact factor (IF) for WJCC as 1.013; IF without journal self cites: 0.991; Ranking: 120 among 165 journals in medicine, general and internal; and Quartile category: Q3.

# **RESPONSIBLE EDITORS FOR THIS ISSUE**

Production Editor: Yan-Xia Xing; Production Department Director: Yun-Xiaojian Wu; Editorial Office Director: Jin-Lei Wang.

NAME OF JOURNAL	INSTRUCTIONS TO AUTHORS
World Journal of Clinical Cases	https://www.wjgnet.com/bpg/gerinfo/204
ISSN	GUIDELINES FOR ETHICS DOCUMENTS
ISSN 2307-8960 (online)	https://www.wjgnet.com/bpg/GerInfo/287
LAUNCH DATE	GUIDELINES FOR NON-NATIVE SPEAKERS OF ENGLISH
April 16, 2013	https://www.wjgnet.com/bpg/gerinfo/240
FREQUENCY	PUBLICATION ETHICS
Semimonthly	https://www.wjgnet.com/bpg/GerInfo/288
EDITORS-IN-CHIEF	PUBLICATION MISCONDUCT
Dennis A Bloomfield, Sandro Vento, Bao-Gan Peng	https://www.wjgnet.com/bpg/gerinfo/208
EDITORIAL BOARD MEMBERS	ARTICLE PROCESSING CHARGE
https://www.wjgnet.com/2307-8960/editorialboard.htm	https://www.wjgnet.com/bpg/gerinfo/242
PUBLICATION DATE	STEPS FOR SUBMITTING MANUSCRIPTS
October 6, 2020	https://www.wjgnet.com/bpg/GerInfo/239
COPYRIGHT	ONLINE SUBMISSION
© 2020 Baishideng Publishing Group Inc	https://www.f6publishing.com

© 2020 Baishideng Publishing Group Inc. All rights reserved. 7041 Koll Center Parkway, Suite 160, Pleasanton, CA 94566, USA E-mail: bpgoffice@wjgnet.com https://www.wjgnet.com



W J C C World Journal of Clinical Cases

Submit a Manuscript: https://www.f6publishing.com

World J Clin Cases 2020 October 6; 8(19): 4286-4302

DOI: 10.12998/wjcc.v8.i19.4286

ISSN 2307-8960 (online)

MINIREVIEWS

# Review of simulation model for education of point-of-care ultrasound using easy-to-make tools

Kyu Chul Shin, Young Rock Ha, Seong-Joon Lee, Jung Hwan Ahn

ORCID number: Kyu Chul Shin 0000-0002-8918-7750; Young Rock Ha 0000-0002-4889-6550; Seong-Joon Lee 0000-0001-9735-6369; Jung Hwan Ahn 0000-0002-4676-1716.

Author contributions: Shin KC and Ahn JH wrote the paper and conceptualization; Shin KC, Ha YR, Lee SJ and Ahn JH collected the data; Ha YR and Lee SJ performed the critical review.

# Conflict-of-interest statement:

Authors declare no conflict of interests for this article.

Open-Access: This article is an open-access article that was selected by an in-house editor and fully peer-reviewed by external reviewers. It is distributed in accordance with the Creative Commons Attribution NonCommercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: htt p://creativecommons.org/licenses /by-nc/4.0/

Manuscript source: Invited manuscript

Received: April 24, 2020 Peer-review started: April 24, 2020 Kyu Chul Shin, Jung Hwan Ahn, Department of Emergency Medicine, Ajou University School of Medicine, Suwon 16499, Gyeonggi-do, South Korea

Young Rock Ha, Department of Emergency Medicine, Bundang Jesaeng Hospital, Seongnam-si 13590, South Korea

Seong-Joon Lee, Department of Neurology, Ajou University School of Medicine, Suwon 16499, Gyeonggi-do, South Korea

Corresponding author: Jung Hwan Ahn, MD, Associate Professor, Attending Doctor, Department of Emergency Medicine, Ajou University School of Medicine, 164, World cup-ro, Yeongtong-gu, Suwon 16499, Gyeonggi-do, South Korea. erdrajh@naver.com

# Abstract

Point-of-care ultrasound (POCUS) is a powerful diagnostic tool and provides treatment guidelines in acute critical settings. However, the limitation of using POCUS is operator dependent. Appropriate and validated training for acquiring and using skills in practice must be conducted before using POCUS in clinical settings in order to keep patients safe. Simulation education models have been introduced as a way to solve and overcome these concerns. However, the commercial simulator with sufficiently secured fidelity is expensive and not always available. This review focused on the inexpensive and easily made simulators for education on POCUS in critical specific situations related to the airway, breathing, circulation, and disability. We introduced the simulators that used non-infectious materials, with easily transportable features, and that had a sonographic appearance reproducibility similar to human tissue. We also introduced the recipe of each simulator in two parts: Materials surrounding disease simulators (surrounding materials) and specific disease simulators themselves (target simulators). This review article covered the following: endotracheal or oesophageal intubation, lung (A-lines, B-lines, lung sliding, and pleural effusions such as hemothorax), central vein access, pericardial fluid (cardiac tamponade), the structure related to the eyes, soft tissue abscess, nerve (regional nerve block), and skull fracture simulators.

Key Words: Critical care; Education; Emergency medicine; Simulation training; Point-ofcare; Ultrasonography

©The Author(s) 2020. Published by Baishideng Publishing Group Inc. All rights reserved.



First decision: August 23, 2020 Revised: September 3, 2020 Accepted: September 12, 2020 Article in press: September 12, 2020 Published online: October 6, 2020

P-Reviewer: El-Shabrawi MH, Ooi L, Tokumine J S-Editor: Zhang L L-Editor: A P-Editor: Xing YX



Core Tip: Inexpensive and easily made simulators for education on point-of-care ultrasound provide a sonographic appearance and resistance similar to human tissue. There were various recipes for making simulators; however, the materials used were similar. These materials were readily available, and the preparation methods were simple. We found that the gelatine with Metamucil and polyvinyl chloride mixtures appear to be most similar to actual human tissue in terms of resistance and sonographic appearance. We introduced each recipe of the simulators in two parts: materials surrounding simulators (surrounding materials) and the specific disease simulators themselves (target simulators).

Citation: Shin KC, Ha YR, Lee SJ, Ahn JH. Review of simulation model for education of pointof-care ultrasound using easy-to-make tools. World J Clin Cases 2020; 8(19): 4286-4302 URL: https://www.wjgnet.com/2307-8960/full/v8/i19/4286.htm DOI: https://dx.doi.org/10.12998/wjcc.v8.i19.4286

# INTRODUCTION

Point-of-care ultrasound (POCUS), which has been implemented as a complementary tool in clinical examinations in the medical field, provides improvements in diagnostic abilities, an increase of confidence, in timely therapeutic decision-making, and in reduction of the complications from invasive procedures<sup>[14]</sup>. It does so particularly in critically unstable patients who need a timely diagnosis and a quick emergent treatment decision. These include cases of hemothorax, pneumothorax, cardiac tamponade, hemoperitoneum, rupture of abdominal aortic aneurysm or dissection, or resuscitation due to shock in an intensive care unit or emergency room. The importance of POCUS in a clinical practice context is well-established<sup>[3,5-7]</sup>.

The most important limitation of ultrasound is known to be operator dependent<sup>[4,5]</sup>. Appropriate and validated training for acquiring skills and transferring these into practice must be conducted before using POCUS in clinical settings, in order to keep patients safe<sup>[3,5,6,8]</sup>. In critical situations in particular, misinterpretations of POCUS findings may threaten patients' safety due to missed or inadequate diagnoses and delays of proper treatment<sup>[3,4]</sup>. Therefore, training on POCUS is a necessary condition before implementing POCUS in critical situations. However, in-field teaching of POCUS and live models with pathology are restricted due to there being few opportunities to use POCUS<sup>[2,9-12]</sup>. There are rarely cases of cardiac tamponade in real clinical practice, and there is also a fear of incurring harm in critical patients, as well as ethical issues and patient safety concerns<sup>[2,9-12]</sup>.

Simulation education models using simulators have presented a way to solve and overcome these concerns, because they provide a safe, stress free, and supportive environment, with excellent reproducibility and fidelity<sup>[2,9-12]</sup>. Furthermore, simulation education models possibly lower the learning curve of skill acquirement<sup>[2,9-12]</sup>. Nowadays, a large number of simulators have been developed and introduced<sup>[2,9,12,13]</sup>. However, the commercial simulator that has a sufficiently secured fidelity is expensive and not always available to trainees<sup>[2,9,12,13]</sup>. Pioneers in the medical field have developed inexpensive and easily made simulators for POCUS educational purposes, and have reported admirable educational achievements<sup>[2,13-16]</sup>. An ideal simulator should be inexpensive, easily made, not be time-consuming to produce, has no infection issues, easily transportable, and has a reproducibility of the sonographic appearance and resistance that is similar to human tissue<sup>[2,9,17]</sup>.

The focus of this review article was to investigate easily made simulators close to or according to the conditions of ideal simulators. The scope of reviewing simulators was related to POCUS use in emergency medicine and critical care medicine, where clinicians must be trained before implementing POCUS in the clinical field to ensure patients' safety. To meet these points of view, this review article is largely divided into two themes. The first part focused on the materials surrounding target simulators such as veins, nerves, or the trachea. The materials surrounding target simulators that have a sonographic appearance and resistance similar to human tissue are also important. This ability to create a surrounding material that is similar to human tissue allows us to make various simulators that are suitable for specific situations. The second part focused on specific target simulator structures related to airway, breathing, circulation, disability, and other situations such as fracture, soft tissue abscesses, and regional nerve blocks. Therefore, this review article covered the following simulators;



endotracheal or oesophageal intubation, lung (A-lines, B-lines, lung sliding, pleural effusion such as hemothorax), central vein access, pericardial fluid (cardiac tamponade), the structure related to eye, soft tissue abscess, nerve (regional nerve block), and skull fracture simulators.

Because of space constraints, the article will not provide the validation methods, results, figures, or creation methods; these are well presented in the original reports on the subject. The review will be limited to the introduction of the recipes or substances used to make the POCUS simulators and the appearance of the ultrasounds. The authors hope that this review article can guide user's searches for specific, easy-tomake made simulators. The authors would like to introduce a website that has a collection of articles and Youtube addresses describing the simple, homemade simulators https://www.ultrasoundtraining.com.au/foamus/diy-phantomcompendium accessed December 27, 2019).

# SIMULATORS SURROUNDING TARGET, SONOGRAPHIC APPEARANCE SIMILAR TO HUMAN TISSUE

An ideal simulator should represent the texture and resistance of human tissue and have a sonographic appearance similar to human tissue<sup>[2,18,19]</sup>. Researchers have been searching for materials that look similar to real human tissue, and using them to come up with simulators<sup>[18-21]</sup>. We would like to introduce some materials used for making simulators to soft human tissue. The method of producing each material cannot be fully described due to space constraints; these are well presented in the original papers on the subjects. The materials, material ratios, methods, and times to produce each simulator are presented briefly in Tables 1 and 2. Table 3 was summarized and compared the pros and cons of various simulation materials.

# Main materials used in making simulators

Gelatine: Gelatine is a common material for ultrasound simulators because it is easy to get, cheap, and has a similar echogenicity to human tissue<sup>[2]</sup>. It is reusable, but has relatively short shelf life and has problems with durability<sup>[22-24]</sup>. Gelatine also requires special ingredients and recipes to make simulators that are more similar to a real tissue and that do not melt at room temperature<sup>[22,24]</sup>. Generally, 100 to 500 g of gelatine is used per 1000 mL of water, and if the proportions of the other additives are high, a small amount of gelatine is used<sup>[13,15,18,22,24-31]</sup>.

Agar: Agar is a promising material for ultrasound simulators because it is not easy to melt or decompose by itself, and does not need to be refrigerated in the process of making it<sup>[32,33]</sup>. However, it has weaknesses in terms of accessibility and cost. Earle et al<sup>[32]</sup> compared several simulators models that used 2.5%-10% agar, and found out that the 5% model had the ideal texture and distensibility. They mixed 750 mL of cold water with 38 g (900 g/cm<sup>2</sup>) of agar gel and stirred the mixture periodically as it boiled<sup>[32]</sup>. They then sprinkled 1 teaspoon of flour into the mixture and stirred it until it homogenized<sup>[32]</sup>. About half of the hot mixture was poured into a container and cooled for 20 min at room temperature<sup>[32]</sup>. Other studies using agar produced simulators using mixtures of similar concentration<sup>[33]</sup>.

Ballistics gel: Ballistics gel was originally developed to replace human tissue in firearms testing<sup>[34,35]</sup>. The density and texture of ballistics gel is similar to human tissues, thus, it is good for ultrasound simulations<sup>[35]</sup>. In practice, its ultrasound images appear very similar to human soft tissue<sup>[34,35]</sup>. It has advantages in many aspects, including durability, cost, and reusable, risk of transmission, and leaking, compared to other materials used in previous studies<sup>[35]</sup>. It can be used repeatedly it is self-sealing, but it becomes slightly yellowish when reused over time<sup>[36]</sup>. The 10% ballistic gelatine (Clear Ballistics, LLC, Fort Smith, AR, United States) which the authors used in their study has a melting point of about 93.3 °C<sup>[34,35]</sup>. Before heating, the ballistic gel was divided into small pieces to decrease the melting time<sup>[34-36]</sup>. Since each manufacturer has a different melting point, it is important to check it well and slowly heat it up and melt it<sup>[34-36]</sup>. Then pour it into a container after melting and let it cool for about 30 to 40 min to harden at room temperature<sup>[34-36]</sup>.

# Special ingredients added to simulators

Additives to increase echogenicity and opacity: To obtain an ultrasound image similar to real human tissue, surrounding materials are made by adding other



Surrounding materials	Ref.	Brand of main materials (A)	Ingredients for echogenicity (B)	Other ingredients (C)	Ratio or amount (A: B: C)	Target simulator
Gelatin based	Lo <i>et al</i> <sup>[18]</sup>	Jell-O gelatin cherry flavor	Sugar-free Metamucil (P and G, Cincinnati, OH, USA) brand psyllium hydrophilic mucilloid fiber	NA	12:4 tablespoons in 1 L water	Soft tissue abscess
	Wilson <i>et al</i> <sup>[53]</sup>	(Kraft Foods, Northfield, IL, United States )				
	Chao et al <sup>[23]</sup>	Unflavored gelatin powder, not mentioned brand name	sugar-free psyllium fiber (Metamucil Sugar Free; P and G)	NA	40 g: 20 g in 250 mL water	Soft tissue abscess, central vein access, foreign body
	Seguin and Tessaro <sup>[29]</sup>	Beef gelatin powder, not mentioned brand name	Orange colored psyllium fiber (Metamucil sugar free; P and G)	NA	90 mL: 60 mL in 500 mL water	Intubation
	Merali <i>et al</i> <sup>[13]</sup>					
	Do and Lee <sup>[28]</sup>	Gelatin, not mentioned brand name	Agar, no mentioned brand name	NA	20 g:1 g in 60 mL water	Lung, A-lines, B-lines, pleural effusion
	Zerth <i>et al</i> <sup>[26]</sup>	Plain gelatin, not mentioned brand name	NA	Green and red food coloring, surgical iodine solution	16 oz (473 mL): several drop of coloring agent + 120 mL iodine in 4 quarts (3785 mL) water	Pericardial effusion
	Cuévas Gonzales <sup>[25]</sup>	Unflavored gelatin, not mentioned brand name	NA	70% ethanol	100 g: 30 mL in 500 mL water	Eye (normal, foreign body, ONS, dislocation of lens, vitreous hemorrhage, retinal detachment)
	Murphy <i>et al</i> <sup>[15]</sup>	Unflavored gelatin, not mentioned brand name	sugar-free psyllium powder (e.g., Metamucil Sugar Free Dietary Fiber Supplement; P and G)	NA	Not mentioned gelatin amount, one tablespoon Metamucil in 250 mL water	ONS
	Soucy <i>et al</i> <sup>[27]</sup>	Sugar-free gelatin, not mentioned brand name	Metamucil (P and G)		20 g: 10 g (1 tablespoon) in 250 mL water	Skull fracture
	Amato <i>et al</i> <sup>[24]</sup>	Unflavored gelatin, not mentioned brand name	NA	NA		Vascular access
	Cheruparambath et al <sup>[22]</sup>	Plain gelatin, not mentioned brand name	Coarsely ground finger millet (Eleusine coracana) flour	NA	15 g: 20 g in 1.2 L water	Vascular access
	Chao et al <sup>[23]</sup>	Unflavored gelatin, not mentioned brand name	Sugar-free psyllium fiber (Metamucil Sugar Free; P and G)	NA	40 g: 20 g in 250 mL water	Vascular access
	Rathbun <i>et al</i> <sup>[39]</sup>	unflavored Knox gelatin (Kraft Foods)	Sugar-free psyllium hydrophilic mucilloid fiber (sugar-free Metamucil; P and G)	NA	4 envelopes of Knox gelatin add 1 tablespoon of Metamucil (2 cups a of cold water)	Nerve block
Agar based	Earle <i>et al</i> <sup>[32]</sup>	Agar, not mentioned brand name		Flour (1 teaspoon) Red dye ethanol	38 g: 750 mL water (for 5% model) Compare 4 different concentrations (10%, 7.5%, 5%, 2.5% by weight)	Vascular access
	Nikitichev <i>et al</i> <sup>[33]</sup>	Agar (A7002; Sigma-Aldrich, St Louis, MO, United States)	NA	NA	5.5% by weight	Vascular access

# Table 1 Ingredients and classification of simulator according to surrounding materials used in simulators

Shin KC et al. Simulation model using easy-to-make tools

Ballistics gel based	Morrow <i>et al</i> <sup>[35]</sup>	10% ballistic gelatin (Clear Ballistics LLC, Fort Smith, AR, United States)	Not required. Main material provided the sonographic appearance similar to human tissue	Clear gel dye if needed	NA	Soft tissue abscess
	Seguin and Tessaro <sup>[29]</sup>	ballistic medical gelatin #3 (Clear Ballistics)	psyllium fiber (Metamucil sugar free; P and G),	ballistic gelatin dye (Clear Ballistics)	250 mL: 5 mL: 45 mL	Intubation
	Young and Kuntz <sup>[38]</sup>	Ballistics gel, plastic Halloween skeleton thorax, plastic, manikin dress-form torso	NA	flesh-colored dye	NA	Pericardial effusion
	Morrow and Broder <sup>[34]</sup>	10% ballistic gelatin (Clear Ballistics LLC)	NA	NA	NA	Vascular access
	Amini et al <sup>[36]</sup>	Clear ballistic gel, not mentioned brand name	NA	NA	NA	Vascular access
Gel and wax based	Campo Dell'orto <i>et al</i> <sup>[20]</sup>	Gel wax (mixed store, 74532 Ilshofen)	Gel (Sonosid 1 L, Asid Bonz GmbH, 71083 Herrenberg, PZN 5362311)	Silicon skin/Thera-Band (Schmidt Sports PHYSIO TAPE, 42699 Solingen, Art.Nr. 111202)	NA	Pericardial effusion
	Daly et al <sup>[37]</sup>	Gel wax, flour, plastic artificial rib cage (axial skeleton)	NA	Silicon skin/Thera-Band (Schmidt Sports PHYSIO TAPE, 42699 Solingen, Art.Nr. 111202)	One gallon (4546 mL) gel wax + three tablespoons of flour	Pericardial effusion

NA: Not applicable; ONS: Optic nerve sheath.

substances<sup>[13,15,18,20,22-30,32-39]</sup>. In gelatine-based simulators, one study used coarsely ground finger millet (Eleusine coracana) flour to increase echogenicity<sup>[24]</sup>. Metamucil was added as an ingredient for echogenicity, however, there was no information of the amount of Metamucil added<sup>[13,15,18,23,27,29,30,39]</sup>. Lo et al<sup>[18]</sup> made a simulators of a skin abscess simulation with gelatine and sugar-free Metamucil. They made several simulators using these two materials in various concentrations and found out a ratio of 12:4 tablespoons in 1 L water of Jell-O brand gelatine to sugar-free Metamucil showed the best results<sup>[18]</sup>. They evaluated the variety of simulators concentrations on both firmness and ultrasound appearance<sup>[18]</sup>. The ratio of 12:4 tablespoons had the best results in both of these aspects<sup>[18]</sup>. There is no validation research on this, as this was the only paper that studied the concentration of these two compounds. Li *et al*<sup>[40]</sup> proposed the use of modified polyvinyl chloride for ultrasound simulators for the first time. They made a novel ultrasound simulator from a mixture of polyvinyl chloride (PVC), mineral oil, diethyl hexyl adipate plasticizer softener, and chalk<sup>[40]</sup>. Pepley et al<sup>[19]</sup> made the modified PVC polymer mixture using three different concentrations for an ultrasound guided vascular access simulators<sup>[19]</sup>. They measured the axial needle force of a needle insertion into the simulators and designed a survey to compare the ultrasound images to real images<sup>[19]</sup>. A total of nine simulators materials were tested including three different concentrations of modified PVC simulators, two kinds of commercial simulators, gelatine, agar, ballistics gel, and cadavers<sup>[19]</sup>. It was found that

Surrounding Material for target						
Surrounding materials	Ref.	Heating	Refrigeration	Time to produce	Target simulator	Material for target simulator
Gelatin based	Lo <i>et al</i> <sup>[18]</sup>	Required	At least 4 h to firm	at least 5 h	Soft tissue abscess	Ultrasound gel
	Wilson <i>et al</i> <sup>[53]</sup>					
	Chao et al <sup>[23]</sup>	Required	NA	NA	Soft tissue abscess, central vein access, foreign body	Pipette bulbs with water, latex rubber tube, NA
	Seguin and Tessaro <sup>[29]</sup>	Required	Required, overnight at 4°C	NA	Intubation. Trachea, esophagus	Holes, the block that cam out when making the hol
	Merali <i>et al</i> <sup>[13]</sup>					
	Do and Lee <sup>[28]</sup>	NA	NA	NA	Lung, A-lines, B-lines, pleural effusion	Dressing material (Medifoam <sup>®</sup> , Mundipharma, South Korea), balloon with wate
	Zerth <i>et al</i> <sup>[26]</sup>	Required	Required for at least 2 h	NA	Pericardial effusion	Balloon with red food coloring water
	Cuévas Gonzales <sup>[25]</sup>	Required	Required for 2 h	NA	Eye (normal, foreign body, ONS, dislocation of lens, vitreous hemorrhage, retinal detachment)	Printing was required using aluminum paper mold on gelatin
	Murphy <i>et al</i> <sup>[15]</sup>	Required	Required	At least 30 min	ONS	Ping-pong ball, Vinyl tub + a stylet of 18-guage needle with metal tip
	Soucy <i>et al</i> <sup>[27]</sup>	Required	Required	NA	Skull fracture	Coconut
	Amato et al <sup>[24]</sup>	Required	Required	NA	Vascular access	Removed plastic tube
	Cheruparambath et al <sup>[22]</sup>	Required	Required	NA	Vascular access	Plastic tube
	Chao et al <sup>[23]</sup>	Required	Required	NA	Vascular access	Latex rubber tube
	Rathbun <i>et al</i> <sup>[39]</sup>	Required	Required	Approximately 2-3 h	Nerve block	Hot dog and straw
Agar based	Earle <i>et al</i> <sup>[32]</sup>	Required	Cooling at room temperature	NA	Vascular access	Penrose drain or long latex balloon
	Nikitichev et al <sup>[33]</sup>	Required	Cooling at room temperature	NA	Vascular access	Removed rod
Ballistics gel based	Morrow <i>et al</i> <sup>[35]</sup>	Required	Cooling at room temperature	Various time according to heating methods	Soft tissue abscess	Water
	Seguin and Tessaro <sup>[29]</sup>	Required	Cooling at room temperature	At least 1 h	Intubation	Marker embedded in surrounding materials
	Young and Kuntz <sup>[38]</sup>	Required	Cooling at room temperature, ice underneath the plastic manikin torso	NA	Pericardial effusion	Ping-pong ball in 250 cc normal saline bag
	Morrow and Broder <sup>[34]</sup>	Required	Cooling at room temperature	Less than 1 h	Vascular access	Latex tube
	Amini et al <sup>[36]</sup>	Required	Cooling at room temperature	NA	Vascular access	Removed plumbing tube
Gel and wax based	Campo Dell'orto <i>et al</i> <sup>[20]</sup>	Required	Not required	1 h	Pericardial effusion	Celluloid table tennis bal with red colored water in balloon
	Daly et al <sup>[37]</sup>	Required	Cooling at room temperature for 1 h	NA	Pericardial effusion	Ping-pong ball in 250 cc normal saline bag

NA: Not applicable; ONS: Optic nerve sheath.

Baisbideng® WJCC | https://www.wjgnet.com

the PVC mixture with a ratio of 11:0:1 (PVC polymer to softener to mineral oil) and ballistics gel provided the best results in terms of needle resistance most similar to cadaveric tissue<sup>[19]</sup>. The survey results showed that a PVC mixture with a 9:1:2 ratio (PVC polymer to softener to mineral oil) was the best performing material with the most realistic image<sup>[19]</sup>. Ballistics gel and agar showed the worst performing materials in the image-comparing survey<sup>[19]</sup>. On both aspects of resistance and ultrasound appearance of the simulators, in addition to the gelatine and Metamucil mixture, we think the PVC mixture model will be able to get the spotlight as a high-performance material. In contrast, agar showed the worst performance as a main material of the simulators, and had a less realistic image and poor needle forces<sup>[41]</sup>.

Surface resistance: In a simulator made of gelatine, the actual skin resistance can be different and it can easily be broken<sup>[20,23]</sup>. Hydrocolloid dressing materials (DuoDERMR CGF® dressings; ConvaTec, Skillman, NJ, United States) were put on the surface of the simulators to avoid disruption of the surrounding material<sup>[23]</sup>. This is a sonolucent material and can be tightly adhered to the surface of surrounding material<sup>[23]</sup>. It has many advantages such as protection from bacteria, avoiding disruption from transducer movement, and can be punctured several times because it increases durability<sup>[23]</sup>. In a pericardial effusion simulation study, Campo Dell'orto et al<sup>[20]</sup> and Daly et al<sup>[37]</sup> covered home-made simulators in silicon skin such as Schmidt Sports PHYSIO TAPE (THERABAND, Akron, OH, United States ), which are used for gymnastic exercise, to make the resistance similar to real skin and prove that it can be used as a skin analog<sup>[20,37]</sup>. While securing skin resistance, it also had the effect of increasing realness<sup>[20,37]</sup>.

Colour: Several substances, such as various colours of dye, can be used to colour surrounding materials<sup>[26,32,38]</sup>. There are different advantages of using transparent materials compared to coloured materials. Transparent materials have external visibility, so they allow a novice to compare real anatomical locations and ultrasound images with the naked eve<sup>[24,34]</sup>. In addition, when the beginner fails a procedure during the simulation, they can confirm the problems themselves<sup>[34]</sup>. For the more advanced learners, coloured simulators can offer more realistic and difficult simulations<sup>[34]</sup>. It is important to add colours to make the models more real, but from an educational point of view, using a transparent simulator to learn relative positions of objects and where the ultrasound probe should go is more important. We thought training using coloured simulators would be of great help as the next step.

# SIMULATORS OF SPECIFIC TARGETS

In this section, we describe the simulators of specific situations in terms of the surrounding material for simulating human tissue (part of surrounding material) and the recipe/materials used to create the simulator for specific disease (part of target simulator). If the recipe was described in the reference, the recipe is introduced briefly. The authors of this review article would like to note that this article does not describe the definition or explanations of the sonographic findings, as that is beyond the scope of this article. However, several sonographic findings were explained as a convenience to readers if sonographic signs are unfamiliar to them. The materials used to make the target simulators are introduced in Table 2. The general method of making easy-made simulator is briefly schematized in Figure 1.

# Related to A- airway

POCUS related to airways is used for detection of cricothyroid membranes or trachea rings for cricothyroidotomy or tracheostomy, tracheal or oesophageal placement of endotracheal tubes, and diaphragmatic movement or lung sliding for confirmation of tracheal placement<sup>[2,5-7,21]</sup>. The simulators related to POCUS of lungs will be introduced in the next section, "related to B- Breathing". To the best of our knowledge and PubMed searches, there have not been any articles related to simulators that detect cricothyroid membrane, trachea, or diaphragmatic movements. The only two articles that we could find for simulators used to confirm the placement of endotracheal tubes (trachea vs oesophagus) were introduced by Merali et al<sup>[13]</sup> and Seguin and Tessaro<sup>[29]</sup>. Merali et al<sup>[13]</sup> reported on the educational effects of the simulator made by Seguin and Tessaro<sup>[29]</sup>. Therefore, the two simulators made by Seguin and Tessaro<sup>[29]</sup> are introduced in this section.

Surrounding material – neck tissue: One simulator was made using beef gelatine



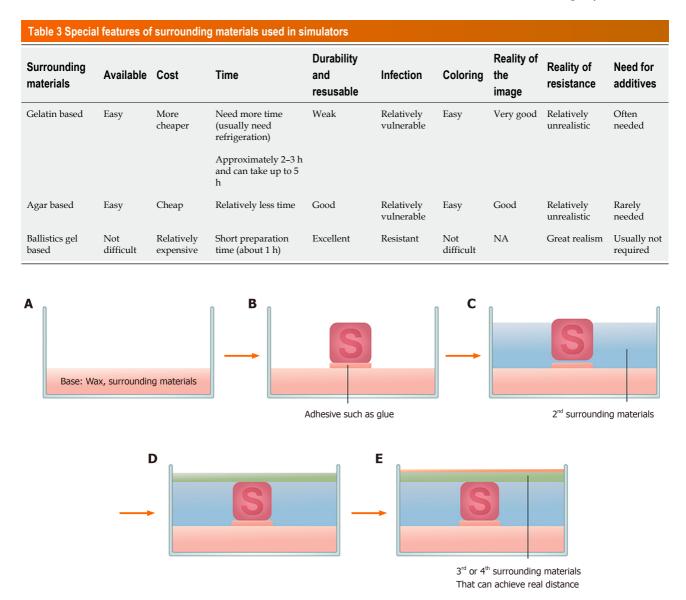


Figure 1 General method to produce the simulator. The surrounding material and specific simulator are used in a container box. A: The surrounding material or wax is used to form a base in a container box prior to the specific simulator insertion; B: When the base material solidifies, the specific simulator is glued with an adhesive; C: The process of inserting the surrounding material; D and E: The surrounding material can be inserted at once or may be split into two depending on the purpose of the project. At this time, air bubble is removed using a spoon.

> powder (90 mL) and 60 mL of orange coloured psyllium fibre (Metamucil sugar free; P and G, Cincinnati, OH, United States) mixed with 500 mL of boiling water<sup>[13,29]</sup>. The mixture was poured into a 1 L container and refrigerated at 4 °C overnight. The other was made using a 250 mL volume of ballistic medical gelatine 3 (Clear Ballistics), psyllium fibre powder (5 mL), and 45 mL of ballistic gelatine dye (Clear Ballistics)<sup>[29]</sup>.

> Target - trachea, oesophagus, and simulation of the placement of endotracheal tubes: Two holes were made in the surrounding materials to simulate a trachea and oesophagus<sup>[13,29]</sup>. The two holes are separated by a gap between the trachea and oesophagus, and are 5 mm from the wall of the margins of the surrounding materials, in order to achieve a real sonographic appearance<sup>[13,29]</sup>. In the first model using beef gelatine, two holes was made by using a 10 mL syringe with the tip cut off, once the surrounding material was completely ready<sup>[13,29]</sup>. In the second model using ballistic gelatine, the holes were made using two felt markers (Crayola)<sup>[29]</sup>. The two felt markers were planted into the surrounding material mixture before solidifying the surrounding materials<sup>[29]</sup>.

> The methods for confirming the placement of the endotracheal tube using POCUS are the static method (double tracheal sign, enhanced tracheal posterior shadow sign, transcricothyroid membrane bullet sign and lung sliding sign) and the dynamic method (snowstorm sign, and shape change of the transcricothyroid membrane from



triangular to bullet sign)<sup>[41,42]</sup>. Among these findings, the double tracheal sign was easily applicable to clinical practice<sup>[42]</sup>. The double tracheal sign simulated an oesophageal intubation<sup>[41,42]</sup>. The presence of air in both the trachea and the oesophageal lumen produces two reverberation artefacts<sup>[42]</sup>.

To simulate the double tracheal sign and bullet sign, Seguin and Tessaro<sup>[29]</sup> used a block that came out when a hole was made using a syringe or felt marker. Therefore, the block was the same material as the surrounding materials. The block in the first simulator that used beef gelatine was made using a syringe while making holes in surrounding material<sup>[29]</sup>. To make the block in the second simulator, the surrounding material mixture (ballistic gelatine) was poured into a 10 mL syringe with the tip cutoff<sup>[29]</sup>. If the block was not planted in the hole that represented the oesophagus, the double trachea sign was shown in ultrasound<sup>[29]</sup>. On the other hand, the simulation with the block in the hole that represented the oesophagus did not show the double tracheal sign, and represented tracheal intubation<sup>[29]</sup>. However, the simulator by Seguin and Tessaro<sup>[29]</sup> has the limitation of not demonstrating all the POCUS signs mentioned above. The stage of making intubation simulator is briefly schematized in Figure 2A.

# Related to B- breathing

POCUS related to breathing such as dyspnoea is used for detection of hemothorax, pleural effusion, pneumothorax, pulmonary embolism, pneumonia, alveolar/interstitial disease, and thoracentesis by using a lung ultrasound<sup>[3,5-7,43]</sup>. The definition of the lung ultrasound findings was not described in this review article because they are out the scope of this review. The definitions of A-line, B-line, lung sliding, and so on in terms of lung ultrasound findings have been well described by Lichtenstein and Meziere<sup>[44]</sup>. Hand ultrasounds were compared to simulated normal lung or pneumothorax models (presence or absence of lung sliding)[45,46]. A gelatinebased model that simulated a normal lung or pneumothorax, pulmonary oedema, or pneumonia (B-lines), and pleural effusion was also tested<sup>[28]</sup>. Another gelatine-based model that simulated the diaphragm, pleural fluid, and lung parenchyma was used in order to train for thoracentesis<sup>[30]</sup>. Sponge and dressing materials in both the hand and gelatine simulators generated A-lines and B-lines<sup>[16,28,47]</sup>.


**Surrounding material – chest wall and ribs:** Simulations using hands are similar to the sonographic appearance of lungs. This is because the parallel bony structure of metacarpal bone is similar to ribs, the soft tissue between the metacarpal bones is similar to the intercostal space, the interface between air and skin is similar to the lung and pleura interface, and similar to the soft tissue of the anterior chest wall<sup>[16,45]</sup>. The gelatine model, involving 20 g of gelatine, 60 mL of water, and 1 g of agar used as an additive, was used to represent the thoracic wall<sup>[28]</sup>. Tree branches were used to make the simulator represent ribs<sup>[28]</sup>. In another report that simulated pleural effusions, a gelatine cage was used to represent the thoracic wall and ribs; however, the recipe was not described<sup>[30]</sup>.

Target - lung sliding, A-line, B-line, and pleural fluid: The dorsal skin of the hand<sup>[16,45]</sup> or the dry polyurethane dressing foam underneath the dorsal hand and gelatine model<sup>[16,28]</sup> represented the pleura and created A-lines. The dorsal skin or dry polyurethane dressing foam with motion from side to side generated lung sliding/seashore sign with A-lines, while no movement generated pneumothorax (absence of lung sliding/stratosphere sign) with A-lines<sup>[16,28,45]</sup>. The educational effects of these lung simulators have been sufficient to teach novices lung ultrasound techniques<sup>[16,28]</sup>. Furthermore, the materials such as sponges and dressing materials were similar to the structure of lungs; therefore, they were able to simulate lungs, A-lines, and Blines<sup>[16,28,30,47]</sup>. Depending on the dry or wet (soaked in water) condition of the sponge and the dressing materials, they were able to create A-lines (dry) or B-lines (wet)<sup>[16,28]</sup>. Pleural effusion was simulated by using a balloon filled with water to contact the bottom of the simulator<sup>[28]</sup>. In the other report that simulated pleural effusion, a sponge and water in a plastic container was used to represent lung and pleural effusion<sup>[30]</sup>. In this simulator, the gelatine simulators that represented the thoracic wall and rib was placed above some plastic wrap, over the sponge and water, in a plastic container<sup>[30]</sup>. This model was able to allow trainees to perform a thoracentesis<sup>[30]</sup>. However, this article did not provide the recipe<sup>[30]</sup>.

# Related to C- circulation

POCUS related to circulation such as shock and vascular access for resuscitation is used for the detection of pericardial effusion such as cardiac tamponade. This can



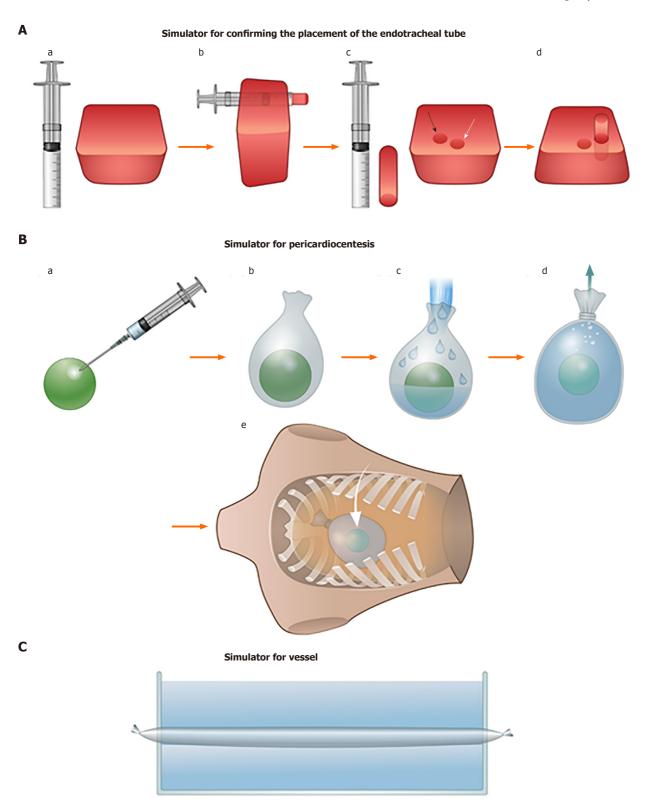

**Figure 2 S specific simulator: The process of making some specific simulator that is complicated to make.** A: Simulator for confirming the placement of the endotracheal tube A (a) beef gelatin powder (90 mL) and 60 mL of orange-colored psyllium fiber (Metamucil sugar free; P and G, Cincinnati, OH, United States) mixed with 500 mL of boiling water in a 1 L container. A 10 mL syringe with the tip cutoff; A (b) The two holes are separated by a gap between the trachea and oesophagus; they are 5 mm from the wall of the margins of the surrounding materials to achieve a real sonographic appearance. Two holes were made using a 10 mL syringe with the tip cutoff once the surrounding material was completely ready; A (c) White arrow shows the trachea. Black arrow shows the oesophagus; A (d) The simulation with the block in the hole that represented the oesophagus did not show the double tracheal sign and showed tracheal intubation. B: Pericardiocentesis simulator B (a) Prick a hole in a ping-pong ball and fill with water. It represents a heart; B (b) The ball then is inserted into a balloon or a 250 cc saline bag, which represents the pericardium; B (c) The balloon is filled with water to show the pericardial fluid; B (d) A spoon or a syringe is used to remove the bubbles; B (e) the simulator is used in an order as in Figure 1, with an artificial rib cage or a dummy. The posterior portion of the artificial rib cage and posterior portion of a dummy are removed so that it can be used as a container. And C: Vessel simulator such as latex or balloon are inserted in process Figure 1.

include pulmonary embolisms, severe dysfunction of left ventricle, regional wall motion abnormalities, evaluation of the aorta (aortic dissection or aortic aneurysm), or determine vascular access for resuscitation by using a vessel ultrasound<sup>[3,5-7,43]</sup>. In terms of focused echocardiography, four articles on cardiac tamponade and pericardiocentesis detection were found<sup>[20,26,37,38]</sup>. In terms of vessel ultrasounds, 11 articles regarding venous catheterization were introduced in this review article<sup>[17,22-24,32-34,36,48-50]</sup>.

Surrounding material - chest wall and soft tissue: In the cardiac tamponade simulator, the surrounding material that represented chest wall was created by Zerth et al<sup>[26]</sup> using gelatine with green food colouring, Campo Dell'sorto et al<sup>[20]</sup> using gel wax, ultrasound gel, silicone skin, three tablespoons of flour, Daly et al[37] using gel wax, flour, artificial rib cage, and latex exercise resistance bands (TheraBand<sup>™</sup>) with a moulding bucket, and Young and Kuntz<sup>[38]</sup> using ballistic gel with flesh-coloured dye and artificial rib cage moulding in a plastic mannequin dress-form torso. What is very interesting is that by using an artificial rib cage or mannequin dress-form torso, the pericardiocentesis simulators were very similar to real thoracic cages<sup>[37,38]</sup>. These simulators were able to provide a more realistic environment by showing external anatomic landmarks, allowing trainees to practice finding the correct positioning of probes and needles<sup>[37,38]</sup>.

Target – heart, pericardial effusion, and pericardial sac: A golf ball<sup>[26]</sup>, water-filled celluloid ball<sup>[20]</sup>, and red-dyed water-filled ping-pong ball<sup>[37,38]</sup> were used to represent heart. The pericardial sac was represented by a 16-inch punching balloon<sup>[26]</sup>, 30 cm balloon<sup>[20]</sup> and 250 mL normal saline bag<sup>[37,38]</sup>. Pericardial effusion in all simulators was represented by red-coloured water<sup>[20,26,37,38]</sup>. To make a thoracic cage with a pericardial sac and heart, it took two to four steps<sup>[20,26,37,38]</sup>. In the first step, a base floor was made for the purpose of a stable base<sup>[20,26,37,38]</sup>. The second step was to fix the correct position of the heart and pericardial sac by pouring a melted gel<sup>[20,26,37,38]</sup> or sonographic gel<sup>[20]</sup> around the heart and pericardial sac while they were centred. The third step was to add the melted gel into a moulding container (Daly *et al*<sup>137</sup>/s method), or mannequin (Young and Kuntz<sup>[38]</sup>'s method) to make the distance between the sternum and the heart simulator similar to a real human. The last step was to cover it with silicone skin to generate a resistance similar to a human and draw the external anatomical marks<sup>[20,37]</sup>. The stage of making pericardiocentesis simulator is briefly schematized in Figure 2B.

Target - vessels: POCUS related to vascular approach techniques has been one of the most widely used practices in emergency departments. Various kinds of materials have been used in several studies to represent vessels; examination glove<sup>[49]</sup>, intravenous infusion sets<sup>[17]</sup>, Penrose drains<sup>[50]</sup>, balloons<sup>[14,32,48]</sup>, silicone tubes<sup>[23,34]</sup>, latex tubes<sup>[23,24]</sup> and plastic tubes<sup>[22]</sup>. Nolting et al<sup>[49]</sup> introduced a simulators made using a water-filled glove finger and SPAM<sup>®</sup>. They used a large straw for driving a tunnel through the entire SPAM<sup>®</sup> block<sup>[49]</sup>. The tunnel played the role of a vessel with a waterfilled glove finger inserted<sup>[49]</sup>. Another method was to plug one end of the tunnel, and fill it with ultrasound gel<sup>[49]</sup>. However, the authors did not discuss which materials were more practical<sup>[49]</sup>. Wells and Goldstein<sup>[17]</sup> suggested using an intravenous infusion setline connected to an IV saline bag as vessel, and polony as the surrounding material. There was no mention of whether it has a similar tactile sensation to real human tissue<sup>[17]</sup>. We thought it seemed easier to replenish a fluid to simulate the "blood." A balloon model was suggested for vascular access<sup>[14,48]</sup>. Fürst *et al*<sup>[14]</sup> used the latex straw balloon to represent the vessel, and chicken breast for the surrounding material. A latex straw balloon, typically used for decoration and sculpturing, was placed into the tunnel (made by Kelly using forceps) in the chicken breast<sup>[14]</sup>. The authors reported that the two combinations had the most similarity to human tissue in echogenicity and compressibility<sup>[14]</sup>. Rippey et al<sup>[48]</sup> also used chicken breast as a material for the simulators and introduced various materials for implementing vessel, including balloons. Various balloon diameters were used to simulate veins; their thin walls simulated the ultrasound appearance and tactile sensation of punctures<sup>[48]</sup>. By adjusting the tension of the balloons, the compressibility made by the ultrasound probe could be determined<sup>[48]</sup>. This can be used to represent a hypovolemic state of a vein by filling it loosely with fluid<sup>[48]</sup>. To express slightly harder feeling of an artery, instead of a balloon, a latex or silicon tube was used<sup>[48]</sup>. The authors made models with varying complexity that could easily be changed, using chicken breasts and various kinds of balloons or silicone tubes<sup>[48]</sup>. They provided video clips where they create various simulators that embody various shapes and types of blood vessels<sup>[48]</sup>.

The materials we introduced earlier were made by inserting "the blood vessel" into the surrounding material. The stage of making vessel simulator is briefly schematized



in Figure 2C. A latex rubber tube was used to melt the gelatine or agar used to make the simulators, instead of using balloons<sup>[23,34]</sup>. We felt that the appearance on ultrasounds, tactile sensation and durability were similar when compared to balloons. Two studies used a Penrose drain, but did not describe its differences compared to other materials<sup>[32,50]</sup>. Sanchez-de-Toledo and Villaverde<sup>[50]</sup> used boneless pieces of chicken breast with a tunnel made using tweezers. A 10 cm Penrose drain was inserted through the tunnel as a vessel<sup>[50]</sup>. Earle et al<sup>[32]</sup> used a Penrose drain and latex balloon as a vessel, with agar as the surrounding material. They recommended that different materials depending on the situation and environment be used to simulate blood vessels<sup>[32]</sup>. Cheruparambath et al<sup>[22]</sup> made a simulator with plastic tubes and gelatine. When compared with other vessel materials such as latex or rubber, using plastic tubes was less realistic and felt more rigid<sup>[22]</sup>. However, it showed advantages in durability and reusability<sup>[22]</sup>.

The materials introduced so far showed more hyperechoic images of vessel walls than actual human tissue does on ultrasounds. Some of the researchers have suggested wall-less simulators in which vessel-mimicking materials are absent<sup>[24,33,36]</sup>. Using these type of simulators, realistic images can be obtained<sup>[24,33,36]</sup>. A simple method to construct wall-less simulators is to remove the rod or tube positioned in the surrounding material, which is composed of gelatine, agar, or ballistic gel<sup>[24,33,36]</sup>. First, they made a hole for the rod on both sides of the plastic chamber and inserted the rod through the chamber<sup>[24,36]</sup>. With the rod in, they then poured the tissue-mimicking material (gelatine, agar, or ballistic gel) and waited until it hardened<sup>[24,36]</sup>. Then they removed the rod through the holes<sup>[24,36]</sup>. However, these methods have been limited in their widespread use due to problems with skill, equipment, and resources<sup>[33]</sup>. In a study proposed by Nikitichev et al<sup>[33]</sup>, a three-dimensional (3D) printer was used to solve these problems and create a simulation model without difficulty. They could also simulate different shapes of blood vessels and various clinical environments without difficulty<sup>[33]</sup>. Using a 3D printer, a rectangular chamber with two holes on each side was designed by two freely available software programs: Blender (Stichting Blender Foundation, Amsterdam, the Netherlands), and FreeCAD (Juergen Riegel, Werner Mayer, and Yorik van Havre; OpenSource, www.freecad.com)<sup>[32]</sup>. Of the two types of 3D printers they used, an additive polymer resin printer (Objet30 Pro; Stratasys, Eden Prairie, MN, United States) using a rigid opaque white or blue material with a gloss finish (VeroWhitePlus RGD835 or VeroBlue; Stratasys) produced better results than the others<sup>[33]</sup>. Agar with a 5.5% weight concentration was used as the surrounding material, and polytetrafluoroethylene (DirectPlastics, Sheffield, England) was used as the rod material<sup>[33]</sup>. The rods had minimal adhesion to the agar<sup>[33]</sup>. We recommend referring to the paper for detailed methods and other details.

# Related to D- disability

POCUS related to disability is used for detection of stroke by transcranial color-coded POCUS via the temporal bone, retinal detachment, a foreign body, vitreous haemorrhage, retrobulbar hematomas, or measurement of the optic nerve sheath diameter (ONSD)<sup>[5,7,25,51,52]</sup>. To the best of our knowledge after several PubMed searches, there have not been any articles related to simulators for transcranial color-coded POCUS. The three articles we did find are related to simulators for ocular POCUS. Jafri *et al*<sup>[51]</sup> made simulators that represented normal anatomy, foreign bodies, vitreous haemorrhages, retinal detachment, and ONSD measurements. The simulator also simulated lens dislocations, and Jafri *et al*<sup>[51]</sup> simulated retrobulbar hematomas. The appearance of the ocular simulator for ONSD measurements made by Murphy et al<sup>[15]</sup> might be closer to the sonographic images of humans than the previous simulators. However, the simulator by Murphy *et al*<sup>[15]</sup> did not represent the anterior structure of the eye (such as lens), and Cuévas Gonzales' method was made using simple methods and simulates more pathologies<sup>[25]</sup>.

Surrounding material soft tissue and vitreous humour: The simulators for ocular POCUS were more focused on the ocular simulators themselves rather than the surrounding materials<sup>[15,25,51]</sup>. Jafri et al<sup>[51]</sup> made eyeball simulators without surrounding materials. Murphy et al<sup>[15]</sup> fixed the ocular simulator using three packs of unflavoured gelatine, one tablespoon of sugar-free psyllium powder (e.g., Metamucil Sugar Free Dietary Fibre Supplement, P and G), and 250 mL water in 473.176 mL plastic cups ( e.g., Solo Plastic Party Cup, Dart Container Corporation, Mason, OH, United States). The surrounding material used in this model<sup>[15]</sup> looked more similar to human tissue than other simulators<sup>[25]</sup>. Cuévas Gonzales<sup>[25]</sup> made a square-shaped block and printed the shape of an eyeball, lens, and iris using aluminium paper. This block was made using 500 mL water, 100 mg of gelatine, and 30 mL of 70% ethanol into a bowl or



Tupperware container<sup>[25]</sup>.

Target - eyeball, retinal detachment, foreign body, vitreous haemorrhage, retrobulbar hematoma, and optic nerve sheath: Using an aluminium paper mould to represent a normal eyeball. They made the eyeball mould, iris mould (anterior chamber), and posterior lens mould using aluminium, and these moulds were printed on the hard gelatine already prepared as described above<sup>[25]</sup>. In addition, a piece of aluminium paper represented a foreign body embedded somewhere in the rear chamber<sup>[25]</sup>. For vitreous haemorrhage, one or two grooves (1 cm to 1.5 cm irregular shapes, using a scalpel) filled with a mix of instant coffee and water were made somewhere in the rear chamber<sup>[25]</sup>. For retinal detachment, they made the retinal detached mould and printed it on the hard gelatine. For lens dislocation, they printed lenses not parallel to irises<sup>[25]</sup>.

Jafri *et al*<sup>[51]</sup> made an eyeball simulator using a round plastic mould similar to the size of the eye. They poured the gelatine mixture and psyllium powder into the round plastic mould<sup>[51]</sup>. The lens appearance was made by a circumferential incision at a 15° angle with the tip of the blade facing the uppermost centre of the eyeball simulator<sup>[51]</sup>. A metal electrocardiograph lead was used to represent the optic nerve sheath (ONS), placed on the bottom of the eyeball simulator<sup>[51]</sup>. For retinal detachment, an incision using a size 15 scalpel was made in the posterior portion of the eyeball simulator<sup>[51]</sup>. For vitreous haemorrhage, a pocket was made by using an 18-gauge needle in the eyeball simulator, and hyperechoic fluid such as glue or casing plaster dissolved in water were injected into the pocket of the eyeball simulator<sup>[51]</sup>. For retrobulbar haemorrhage, the structure filled with fluid was placed under the posterior portion of the eyeball simulator<sup>[51]</sup>. Of course, it was possible to represent ocular foreign bodies by inserting various materials such as plastic or metal materials into the model as well<sup>[51]</sup>.

Murphy et al<sup>[15]</sup> focused on creating an ocular simulator to represent ONS. To represent ONS, a 7 cm clear vinyl tube with a similar width to the optic nerve and a stylet of an 18-guage needle with a metal tip were used<sup>[15]</sup>. Eyeballs were represented by using a 40 mm diameter ping-pong ball filled with water<sup>[15]</sup>. To make the eyeball simulator, the ONS simulator was inserted in a 473.176 mL plastic cup through a round hole in the middle of the bottom (created using a drill), and then the surrounding material mentioned above was poured into the plastic cup at the level of the metal tip<sup>[15]</sup>. After hardening the surrounding materials in the plastic cup, the pingpong ball filled with water was placed on the ONS simulator, and then the surrounding materials were poured into the plastic cup on top of the ping-pong ball<sup>[15]</sup>.

#### Others - fracture, regional nerve block, and soft tissue abscess

Musculoskeletal ultrasounds are useful tools to diagnose skin infections and fractures, as are other imaging modalities such as plain radiography, computed tomography, and magnetic resonance imaging<sup>[3,7]</sup>. Ultrasounds in particular are superior in situations where other modalities are not suitable or accessible<sup>[3,7]</sup>. Physical examination is limited when it comes to distinguishing between soft tissue abscess and cellulitis<sup>[18,53]</sup>. POCUS has been reported to be a powerful diagnosis tool to diagnose soft tissue abscess or cellulitis, as well as help make treatment decisions<sup>[18,49,53]</sup>. POCUS is useful when it comes to regional nerve blocks in emergency departments, including forearm nerve blocks, brachial plexus blocks, serratus anterior plane blocks, and fascia iliaca compartment blocks<sup>[31]</sup>. In this section, we described the simulators related to musculoskeletal ultrasounds, including fractures, regional nerve blocks, and soft tissue abscesses.

**Fracture:** A paediatric skull fracture simulator has been introduced by Soucy *et al*<sup>[27]</sup>. Coconuts were used to represent paediatric skulls<sup>[27]</sup>. A coconut was cut into puzzlelike pieces, similar to suture lines and skull fractures. The motorized saw has been reported to be the best cutting tool<sup>[27]</sup>. The fixation among coconut pieces was done by applying the epoxy gel to the inner surface of the coconut pieces, in order to maintain the echogenic properties of the coconut and the suture line<sup>[27]</sup>. To waterproof the coconut simulator, a waterproof shellac was applied in two thin coats<sup>[27]</sup>. To stabilize the simulator, two 1.5 cm long pieces of round wooden dowels were attached to the convexity of the simulator, 1 cm from the midsagittal suture on each side<sup>[27]</sup>. To make soft tissue that is similar to the human cranium, the coconut simulator was submerged in a round bowl filled with a sugar-free gelatine and Metamucil (P and G) solution several times<sup>[27]</sup>. The completed skull fracture simulator looked like a lollipop, and was able to represent sonographic suture lines and fracture lines<sup>[27]</sup>.



Regional nerve blocks: The nerve block simulators were made by various materials to simulate nerves such as metal rods<sup>[54]</sup>, wooden dowels<sup>[55]</sup>, electrical wire<sup>[55]</sup>, shoelaces<sup>[17]</sup>, yarns<sup>[31]</sup>, Computer Aided Translation 5 Ethernet cables<sup>[56]</sup>, spaghetti noodles<sup>[57]</sup>, and straw embedded hot-dogs<sup>[39]</sup>. Sparks et al<sup>[56]</sup> made a nerve block simulator using pressure-injected ultrasound gel with a Computer Aided Translation 5 Ethernet cable representing a nerve. This was a good simulator, which produced the "donut sign" of hypoechoic fluid around the nerve, but due to the metallic properties, posterior shadowing appeared<sup>[56]</sup>. Micheller et al<sup>[57]</sup> used a long tubular balloon (28.5 cm × 0.9 cm in dimension, commonly used for balloon animals) filled with water and five pieces of flour spaghetti noodles presenting the femoral nerves. The regular-thickness spaghetti noodles created the "honeycomb" appearance that is typical in sonographic findings of large nerves<sup>[57]</sup>. The visualized hydrodissection could be confirmed in this simulator, but after three or four blockade attempts, it gradually lost hydrodissection confirmation<sup>[57]</sup>. Naraghi et al<sup>[31]</sup> introduced meat glue (transglutaminase) as a binding component to create fascial planes, in order to represent hydrodissection, by injecting fluid through a needle. To substitute the role of the fascia, the meat glue was sprinkled between a pork tenderloin that was cut into the required muscle shape<sup>[31]</sup>. Chunky yarn soaked in ultrasound gel was used as nerves[31]. Several strands of yarn could be twisted together for representing thicker nerves<sup>[31]</sup>. In addition, thin size straws filled with ultrasound gel that were occluded on both ends with chewing gum were used to represent the interscalene brachial plexus<sup>[31]</sup>. Situ-LaCasse et al<sup>[58]</sup> used a nerve simulator that was created by inserting angel hair pasta into long balloons and filling them with a dash of water to educate medical students in ultrasound-guided peripheral nerve block procedures. Previous research models had the disadvantage of having posterior shadowing, so it was hard to confirm the needle tip behind "the nerve"<sup>[39]</sup>. Rathbun et al<sup>[39]</sup> created a realistic nerve block simulator with minimal posterior acoustic shadowing compared to previous simulators. The authors used a gelatine-Metamucil mixture and a hot dog embedded with a straw in the core<sup>[39]</sup>. The centre of the hot dog with a straw represented a nerve<sup>[39]</sup>. They could simulate various nerve target models by adjusting the depth of the surrounding tissue, done by regulating the amount of the mixture used per layer, and the size of the nerves could be varied using different straw sizes<sup>[39]</sup>. However, it was difficult to confirm the perinerve infiltration or hydrodissection, since most of the injectate proceeded retrogradely along the path of the needle<sup>[39]</sup>.

**Soft tissue abscesses:** The simulators for soft tissue abscess have been made using various liquids such as ketchup<sup>[18]</sup>, lotion<sup>[18]</sup>, mayonnaise<sup>[18]</sup>, candies with liquid centres<sup>[18]</sup>, honey<sup>[18]</sup>, ultrasound gel<sup>[18,49]</sup>, mustard<sup>[18]</sup>, lubricating jelly<sup>[18]</sup>, water<sup>[35]</sup>, tapioca with olive oil<sup>[53]</sup>, vanilla pudding with olive oil<sup>[53]</sup>, and salad dressing<sup>[18]</sup> in cystic like materials such as pipette bulbs sealed using kitchen clips<sup>[23]</sup>, or balloons<sup>[35,53]</sup>. Ultrasound gel has been reported to have the best performance in terms of sonographic appearance and injection into surrounding materials, as well as drainage using syringes<sup>[18,49]</sup>. The surrounding materials used in soft tissue abscess were as follows; SPAM<sup>®[49]</sup>, 10% ballistic gelatine<sup>[35]</sup>, and gelatine with sugar-free Metamucil<sup>[18,23,53]</sup>.

# RECOMMENDATIONS OF THE AUTHORS

Each way of producing the simulator has both advantages and disadvantages. Hence, our recommendation on a single method may not be appropriate. In addition, the purpose and circumstances may be different for each institute. The surrounding material recommendation maybe an exception that is based on the review and the authors' experience, which is the closest in human tissue that has a similar sonographic appearance and resistance. In addition, limited studies have used a specific target simulator. It is difficult to recommend one method since all methods have advantages and are similar to real ones. It is best to choose a specific target simulator depending on the structure they wish to focus on. Using hand ultrasound, which is the easiest approach, a dummy that is very similar to real and a pericardiocentesis model that uses a rib cage and mannequin have been developed.

Based on the current review and our experience, we would like to recommend the below simulation models. It is important to add colors to make the models more real, but from an educational point of view, using a transparent simulator to learn the relative positions of the objects and where the ultrasound probe should be placed is more important. We thought that training using transparent simulators would be



more helpful to novices who are unfamiliar with ultrasound and that colored simulators would be of great help in the subsequent steps. To increase the best ultrasound appearance, we recommend the following additives: the gelatin and Metamucil mixture<sup>[18]</sup> (12:4 tablespoons) and the PVC mixture model<sup>[19]</sup> [9:1:2 ratio (PVC polymer to softener to mineral oil)]. We also recommend the modified PVC polymer model in terms of mimicking the skin resistance. If the PVC mixture ratio is practically complicated, we also recommend using hydrocolloid dressing materials (DuoDERMR CGF® dressings; ConvaTec, Skillman, NJ, United States) to secure the skin resistance or Schmidt Sports PHYSIO TAPE (THERABAND, Akron, OH, United States), which also has the effect of increasing realness<sup>[23]</sup> with the gelatin and Metamucil mixture<sup>[18]</sup> (12:4 tablespoons). Table 3 shows the advantages or disadvantages of the surrounding materials. For specific simulation models, we recommend using hand ultrasound with a dressing material form as a practical way of familiarizing the trainee with the basic sonographic appearance and physiology of the lungs<sup>[16,45]</sup>. For the pericardial model, we recommend the mannequin (Young et al<sup>[38]</sup>'s method), which is a compilation of the previous methods. For mimicking the orbit, different recommendations may be provided depending on the tissue that was focused to visualize. Various types of materials have been used in several studies to represent vessels. It is our opinion that for most cases, veins may be best represented by balloons and arteries may be best represented by silicon or latex pipes. Although innovative, 3D printing models do not seem practical so far. Finally, we wish to point out that the aforementioned models have their own strongpoints, and the recommendations provided are only based on our opinion.

# CONCLUSION

Inexpensive and easily made simulators for education on POCUS provide sonographic appearances and resistance similar to human tissue. There were various recipes for making simulators; however, the materials used were similar. These materials were readily available, and the methods of making the simulators were simple. However, one thing to keep in mind is that injuries such as burns can occur in the process of making simulators, particularly during boiling or melting processes using stoves. The sonographic appearance of the easily made simulators that introduced in this review article might not be completely similar to human tissue. However, since the primary goal in procedural education is target acquisition, coordination between targets and needles, and needle finding, the sonographic appearance of these simulators did not affect the utility of the models in achieving those goals. In terms of resistance and sonographic appearance of the simulators, the gelatine and Metamucil mixtures and polyvinyl chloride mixtures appear to be most similar to the characteristics of actual human tissue. They get the spotlight of being a high-performance material.

There are some limitations in this review article. We did not describe the validation methods used to investigate the educational effects, and we did not insert the figures that introduced the creation processes of the simulators due to space constraints. We hope that this review article provides guidance when searching for specific simulators. If the descriptions of the processes to make the simulators were insufficient, it would be helpful to refer to the figures in the original articles.

# REFERENCES

- Gogalniceanu P, Sheena Y, Kashef E, Purkayastha S, Darzi A, Paraskeva P. Is basic emergency ultrasound training feasible as part of standard undergraduate medical education? J Surg Educ 2010; 67: 152-156 [PMID: 20630425 DOI: 10.1016/j.jsurg.2010.02.008]
- 2 Farjad Sultan S, Shorten G, Iohom G. Simulators for training in ultrasound guided procedures. Med Ultrason 2013; 15: 125-131 [PMID: 23702502 DOI: 10.11152/mu.2013.2066.152.sfs1gs2]
- Gleeson T, Blehar D. Point-of-Care Ultrasound in Trauma. Semin Ultrasound CT MR 2018; 39: 374-383 [PMID: 30070230 DOI: 10.1053/j.sult.2018.03.007]
- Andersen CA, Holden S, Vela J, Rathleff MS, Jensen MB. Point-of-Care Ultrasound in General Practice: A Systematic Review. Ann Fam Med 2019; 17: 61-69 [PMID: 30670398 DOI: 10.1370/afm.2330]
- 5 American College of Emergency Physicians. Emergency ultrasound guidelines. Ann Emerg Med 2009; 53: 550-570 [PMID: 19303521 DOI: 10.1016/j.annemergmed.2008.12.013]
- Mayo PH, Beaulieu Y, Doelken P, Feller-Kopman D, Harrod C, Kaplan A, Oropello J, Vieillard-Baron A, 6 Axler O, Lichtenstein D, Maury E, Slama M, Vignon P. American College of Chest Physicians/La Société de Réanimation de Langue Française statement on competence in critical care ultrasonography. Chest 2009; 135: 1050-1060 [PMID: 19188546 DOI: 10.1378/chest.08-2305]



- Bobbia X, Abou-Badra M, Hansel N, Pes P, Petrovic T, Claret PG, Lefrant JY, de La Coussaye JE; 7 Winfocus France Group. Changes in the availability of bedside ultrasound practice in emergency rooms and prehospital settings in France. Anaesth Crit Care Pain Med 2018; 37: 201-205 [PMID: 28826982 DOI: 0.1016/j.accpm.2017.06.008
- Wong I, Jayatilleke T, Kendall R, Atkinson P. Feasibility of a focused ultrasound training programme for 8 medical undergraduate students. Clin Teach 2011; 8: 3-7 [PMID: 21324064 DOI: 10.1111/j.1743-498X.2010.00416.x]
- 9 Nicholson RA, Crofton M. Training phantom for ultrasound guided biopsy. Br J Radiol 1997; 70: 192-194 [PMID: 9135447 DOI: 10.1259/bjr.70.830.9135447]
- Naik VN, Matsumoto ED, Houston PL, Hamstra SJ, Yeung RY, Mallon JS, Martire TM. Fiberoptic 10 orotracheal intubation on anesthetized patients: do manipulation skills learned on a simple model transfer into the operating room? Anesthesiology 2001; 95: 343-348 [PMID: 11506104 DOI: 10.1097/00000542-200108000-00014
- Gallagher AG, Ritter EM, Champion H, Higgins G, Fried MP, Moses G, Smith CD, Satava RM. Virtual 11 reality simulation for the operating room: proficiency-based training as a paradigm shift in surgical skills training. Ann Surg 2005; 241: 364-372 [PMID: 15650649 DOI: 10.1097/01.sla.0000151982.85062.80]
- Niazi AU, Ramlogan R, Prasad A, Chan VW. A new simulation model for ultrasound-aided regional 12 anesthesia. Reg Anesth Pain Med 2010; 35: 320-321 [PMID: 20460975 DOI: 10.1097/AAP.0b013e3181df226b]
- Merali HS, Tessaro MO, Ali KQ, Morris SK, Soofi SB, Ariff S. A novel training simulator for portable 13 ultrasound identification of incorrect newborn endotracheal tube placement - observational diagnostic accuracy study protocol. BMC Pediatr 2019; 19: 434 [PMID: 31722685 DOI: 10.1186/s12887-019-1717-y]
- Fürst RV, Polimanti AC, Galego SJ, Bicudo MC, Montagna E, Corrêa JA. Ultrasound-Guided Vascular 14 Access Simulator for Medical Training: Proposal of a Simple, Economic and Effective Model. World J Surg 2017; **41**: 681-686 [PMID: 27766394 DOI: 10.1007/s00268-016-3757-x]
- Murphy DL, Oberfoell SH, Trent SA, French AJ, Kim DJ, Richards DB. Validation of a Low-cost Optic 15 Nerve Sheath Ultrasound Phantom: An Educational Tool. J Med Ultrasound 2017; 25: 96-100 [PMID: 30065467 DOI: 10.1016/j.jmu.2017.01.003]
- Park EJ, Yoon YT, Hong CK, Ha YR, Ahn JH. Randomized, noninferiority study between video vs hand 16 ultrasound with wet foam dressing materials to simulate B-lines in lung ultrasound: A CONSORT-compliant article. Medicine (Baltimore) 2017; 96: e7642 [PMID: 28746228 DOI: 10.1097/MD.00000000007642]
- Wells M, Goldstein L. The polony phantom: a cost-effective aid for teaching emergency ultrasound 17 procedures. Int J Emerg Med 2010; 3: 115-118 [PMID: 20606820 DOI: 10.1007/s12245-009-0156-1]
- 18 Lo MD, Ackley SH, Solari P. Homemade ultrasound phantom for teaching identification of superficial soft tissue abscess. Emerg Med J 2012; 29: 738-741 [PMID: 21946182 DOI: 10.1136/emermed-2011-200264]
- Pepley DF, Sonntag CC, Prabhu RS, Yovanoff MA, Han DC, Miller SR, Moore JZ. Building Ultrasound 19 Phantoms With Modified Polyvinyl Chloride: A Comparison of Needle Insertion Forces and Sonographic Appearance With Commercial and Traditional Simulation Materials, Simul Healthc 2018; 13: 149-153 [PMID: 29620705 DOI: 10.1097/SIH.000000000000302]
- 20 Campo Dell'orto M, Hempel D, Starzetz A, Seibel A, Hannemann U, Walcher F, Breitkreutz R. Assessment of a low-cost ultrasound pericardiocentesis model. Emerg Med Int 2013; 2013: 376415 [PMID: 24288616 DOI: 10.1155/2013/376415]
- You-Ten KE, Siddiqui N, Teoh WH, Kristensen MS. Point-of-care ultrasound (POCUS) of the upper 21 airway. Can J Anaesth 2018; 65: 473-484 [PMID: 29349733 DOI: 10.1007/s12630-018-1064-8]
- Cheruparambath V, Sampath S, Deshikar LN, Ismail HM, Bhuvana K. A low-cost reusable phantom for 22 ultrasound-guided subclavian vein cannulation. Indian J Crit Care Med 2012; 16: 163-165 [PMID: 23188960 DOI: 10.4103/0972-5229.1020971
- Chao SL, Chen KC, Lin LW, Wang TL, Chong CF. Ultrasound phantoms made of gelatin covered with 23 hydrocolloid skin dressing. J Emerg Med 2013; 45: 240-243 [PMID: 23399392 DOI: 10.1016/j.jemermed.2012.11.022
- Amato ACM, Freitas SLd, Veloso PM, Correia TCV, Santos RVd, Amato SJdTA. Gelatin model for 24 training ultrasound-guided puncture. Jornal Vascular Brasileiro 2015; 14: 200-204 [DOI: 10.1590/1677-5449.0088]
- Burgueño P, Garmilla P. 8th WINFOCUS World Congress on Ultrasound in Emergency and Critical Care 25 Barcelona, Spain. 20-23 October 2012. Abstracts. Crit Ultrasound J 2012; 4 Suppl 1: A1-A30 [PMID: 23282109 DOI: 10.1186/2036-7902-4-S1-A26
- Zerth H, Harwood R, Tommaso L, Girzadas DV Jr. An inexpensive, easily constructed, reusable task trainer 26 for simulating ultrasound-guided pericardiocentesis. J Emerg Med 2012; 43: 1066-1069 [PMID: 21925818 DOI: 10.1016/i.jemermed.2011.05.0661
- Soucy ZP, Mills L, Rose JS, Kelley K, Ramirez F, Kuppermann N. Creation of a High-fidelity, Low-cost 27 Pediatric Skull Fracture Ultrasound Phantom. J Ultrasound Med 2015; 34: 1473-1478 [PMID: 26206834 DOI: 10.7863/ultra.34.8.1473]
- Do HH, Lee S. A Low-Cost Training Phantom for Lung Ultrasonography. Chest 2016; 150: 1417-1419 28 [PMID: 27938759 DOI: 10.1016/j.chest.2016.09.033]
- 29 Seguin J, Tessaro MO. A Simple, Inexpensive Phantom Model for Intubation Ultrasonography Training. Chest 2017; 151: 1194-1196 [PMID: 28483123 DOI: 10.1016/j.chest.2017.02.014]
- Wilson J, Myers C, Lewiss RE. A low-cost, easy to make ultrasound phantom for training healthcare 30 providers in pleural fluid identification and task simulation in ultrasound-guided thoracentesis, Vis J Emerg Med 2017; 80-81 [DOI: 10.1016/j.visj.2017.03.010]
- Naraghi L, Lin J, Odashima K, Buttar S, Haines L, Dickman E. Ultrasound-guided regional anesthesia 31 simulation: use of meat glue in inexpensive and realistic nerve block models. BMC Med Educ 2019; 19: 145 [PMID: 31092233 DOI: 10.1186/s12909-019-1591-1]
- Earle M, Portu G, DeVos E. Agar ultrasound phantoms for low-cost training without refrigeration. Afr J 32 Emerg Med 2016; 6: 18-23 [PMID: 30456059 DOI: 10.1016/j.afjem.2015.09.003]
- 33 Nikitichev DI, Barburas A, McPherson K, Mari JM, West SJ, Desjardins AE. Construction of 3-Dimensional



Printed Ultrasound Phantoms with Wall-less Vessels. J Ultrasound Med 2016; **35**: 1333-1339 [PMID: 27162278 DOI: 10.7863/ultra.15.06012]

- 34 Morrow DS, Broder J. Cost-effective, Reusable, Leak-resistant Ultrasound-guided Vascular Access Trainer. J Emerg Med 2015; 49: 313-317 [PMID: 26093938 DOI: 10.1016/j.jemermed.2015.04.005]
- 35 Morrow DS, Cupp JA, Broder JS. Versatile, Reusable, and Inexpensive Ultrasound Phantom Procedural Trainers. J Ultrasound Med 2016; 35: 831-841 [PMID: 26969595 DOI: 10.7863/ultra.15.04085]
- 36 Amini R, Kartchner JZ, Stolz LA, Biffar D, Hamilton AJ, Adhikari S. A novel and inexpensive ballistic gel phantom for ultrasound training. *World J Emerg Med* 2015; 6: 225-228 [PMID: 26401186 DOI: 10.5847/wjem.j.1920-8642.2015.03.012]
- 37 Daly R, Planas JH, Edens MA. Adapting Gel Wax into an Ultrasound-Guided Pericardiocentesis Model at Low Cost. West J Emerg Med 2017; 18: 114-116 [PMID: 28116020 DOI: 10.5811/westjem.2016.10.31506]
- 38 Young TP, Kuntz HM. Modification of Daly's Do-it-yourself, Ultrasound-guided Pericardiocentesis Model for Added External Realism. West J Emerg Med 2018; 19: 465-466 [PMID: 29760841 DOI: 10.5811/westjem.2018.2.37736]
- 39 Rathbun KM, Brader WT, Norbury JW. A Simple, Realistic, Inexpensive Nerve Phantom. J Ultrasound Med 2019; 38: 2203-2207 [PMID: 30575067 DOI: 10.1002/jum.14905]
- 40 Li W, Belmont B, Shih A. Design and Manufacture of Polyvinyl Chloride (PVC) Tissue Mimicking Material for Needle Insertion. *Procedia Manufacturing* 2015; 866-878 [DOI: 10.1016/j.promfg.2015.09.078]
- 41 Abbasi S, Farsi D, Zare MA, Hajimohammadi M, Rezai M, Hafezimoghadam P. Direct ultrasound methods: a confirmatory technique for proper endotracheal intubation in the emergency department. *Eur J Emerg Med* 2015; 22: 10-16 [PMID: 24441884 DOI: 10.1097/MEJ.00000000000108]
- 42 Hoffmann B, Gullett JP, Hill HF, Fuller D, Westergaard MC, Hosek WT, Smith JA. Bedside ultrasound of the neck confirms endotracheal tube position in emergency intubations. *Ultraschall Med* 2014; 35: 451-458 [PMID: 25014479 DOI: 10.1055/s-0034-1366014]
- 43 Ha YR, Toh HC. Clinically integrated multi-organ point-of-care ultrasound for undifferentiated respiratory difficulty, chest pain, or shock: a critical analytic review. *J Intensive Care* 2016; 4: 54 [PMID: 27529030 DOI: 10.1186/s40560-016-0172-1]
- 44 Lichtenstein DA, Mezière GA. Relevance of lung ultrasound in the diagnosis of acute respiratory failure: the BLUE protocol. Chest 2008; 134: 117-125 [PMID: 18403664 DOI: 10.1378/chest.07-2800]
- 45 Shokoohi H, Boniface K. Hand ultrasound: a high-fidelity simulation of lung sliding. Acad Emerg Med 2012; 19: E1079-E1083 [PMID: 22978736 DOI: 10.1111/j.1553-2712.2012.01431.x]
- 46 Ahn JH, Jeon J, Toh HC, Noble VE, Kim JS, Kim YS, Do HH, Ha YR. SEARCH 8Es: A novel point of care ultrasound protocol for patients with chest pain, dyspnea or symptomatic hypotension in the emergency department. *PLoS One* 2017; 12: e0174581 [PMID: 28355246 DOI: 10.1371/journal.pone.0174581]
- 47 Blüthgen C, Sanabria S, Frauenfelder T, Klingmüller V, Rominger M. Economical Sponge Phantom for Teaching, Understanding, and Researching A- and B-Line Reverberation Artifacts in Lung Ultrasound. J Ultrasound Med 2017; 36: 2133-2142 [PMID: 28626903 DOI: 10.1002/jum.14266]
- 48 Rippey JC, Blanco P, Carr PJ. An affordable and easily constructed model for training in ultrasound-guided vascular access. J Vasc Access 2015; 16: 422-427 [PMID: 26349885 DOI: 10.5301/jva.5000384]
- 49 Nolting L, Hunt P, Cook T, Douglas B. An Inexpensive and Easy Ultrasound Phantom: A Novel Use for SPAM. J Ultrasound Med 2016; 35: 819-822 [PMID: 26939600 DOI: 10.7863/ultra.14.06023]
- 50 Sanchez-de-Toledo J, Villaverde I. Advanced Low-Cost Ultrasound-Guided Vascular Access Simulation: The Chicken Breast Model. *Pediatr Emerg Care* 2017; 33: e43-e45 [PMID: 26599465 DOI: 10.1097/PEC.00000000000620]
- 51 Jafri F, Runde D, Saul T, Lewiss RE. An inexpensive and easy simulation model of ocular ultrasound that mimics normal anatomy as well as abnormal ophthalmologic conditions. *J Ultrasound Med* 2011; 30: 569-573 [PMID: 21460157 DOI: 10.7863/jum.2011.30.4.569]
- 52 Whitson MR, Mayo PH. Ultrasonography in the emergency department. *Crit Care* 2016; 20: 227 [PMID: 27523885 DOI: 10.1186/s13054-016-1399-x]
- 53 Wilson J, Ng L, Browne V, Lewiss RE. An Easy-to-Make, Low-Cost Ultrasound Phantom for Simulation Training in Abscess Identification and Aspiration. *J Ultrasound Med* 2017; 36: 1241-1244 [PMID: 28258634 DOI: 10.7863/ultra.16.04077]
- 54 van Geffen GJ, Mulder J, Gielen M, van Egmond J, Scheffer GJ, Bruhn J. A needle guidance device compared to free hand technique in an ultrasound-guided interventional task using a phantom. *Anaesthesia* 2008; 63: 986-990 [PMID: 18540929 DOI: 10.1111/j.1365-2044.2008.05524.x]
- 55 Pollard BA. New model for learning ultrasound-guided needle to target localization. *Reg Anesth Pain Med* 2008; 33: 360-362 [PMID: 18675749 DOI: 10.1016/j.rapm.2008.02.006]
- 56 Sparks S, Evans D, Byars D. A low cost, high fidelity nerve block model. Crit Ultrasound J 2014; 6: 12 [PMID: 25411589 DOI: 10.1186/s13089-014-0012-2]
- 57 Micheller D, Chapman MJ, Cover M, Porath JD, Theyyunni N, Kessler R, Huang R. A low-fidelity, high-functionality, inexpensive ultrasound-guided nerve block model. *CJEM* 2017; 19: 58-60 [PMID: 27478036 DOI: 10.1017/cem.2016.335]
- 58 Situ-LaCasse EH, Amini R, Bain V, Acuña J, Samsel K, Weaver C, Valenzuela J, Pratt L, Patanwala AE, Adhikari S. Performance of Ultrasound-guided Peripheral Nerve Blocks by Medical Students After One-day Training Session. *Cureus* 2019; 11: e3911 [PMID: 30931182 DOI: 10.7759/cureus.3911]

Zaishidena® WJCC | https://www.wjgnet.com



# Published by Baishideng Publishing Group Inc 7041 Koll Center Parkway, Suite 160, Pleasanton, CA 94566, USA Telephone: +1-925-3991568 E-mail: bpgoffice@wjgnet.com Help Desk: https://www.f6publishing.com/helpdesk https://www.wjgnet.com

